

Watch Tutorial Video

Scan the QR Code to watch our tutorial on how to use the device for better treatment results.






How-to-Use QR Code

EM50 How To Use

Official Website Email: support.us@jmoon.com
US AMAZON Email: amazonsupport@jmoon.com

Tel: +1 855-681-8599
Mon - Fri: 10AM-7PM EST
www.jmoonglobal.com

About the Device

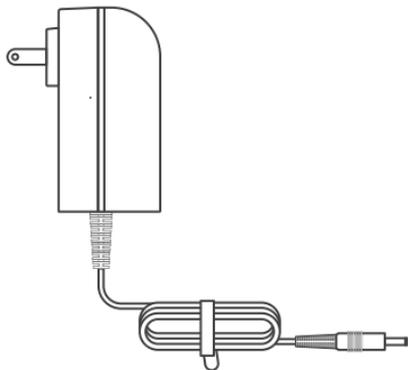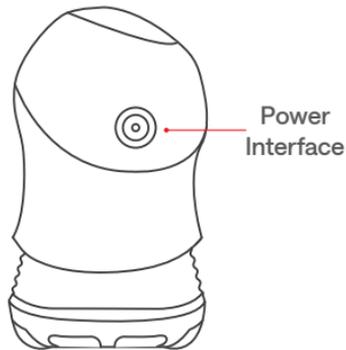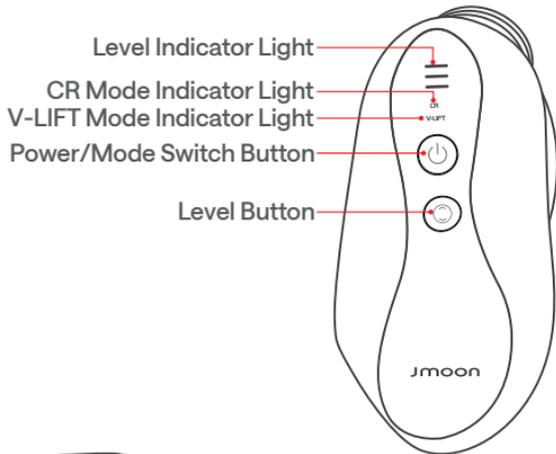


Power Adapter




JMOON Gel



User Manual

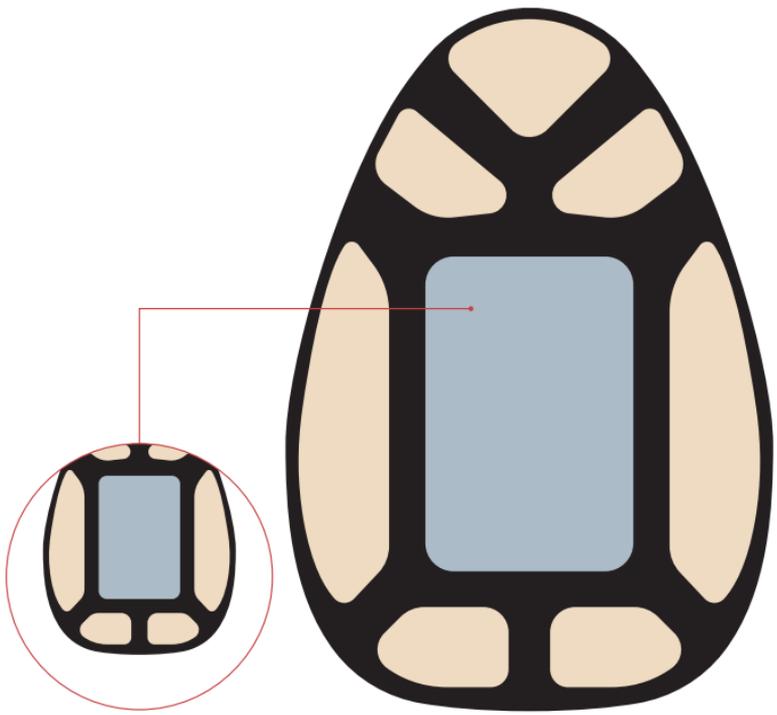


Light Outlet

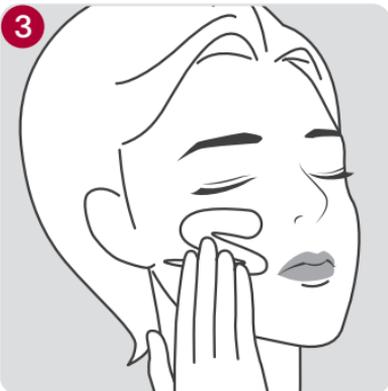


Apply The Gel

- The JMOON Soothing Moisturizing Gel **MUST** be applied before use to achieve your desired results. It is not optional.
 - Apply a **THICK** layer of gel to the operation area **only**.
 - Reapply gel for each operation area.
 - If the device does **NOT** glide smoothly, reapply more gel.
- * Use a **generous amount of gel for a thick layer in order to get optimal results.**

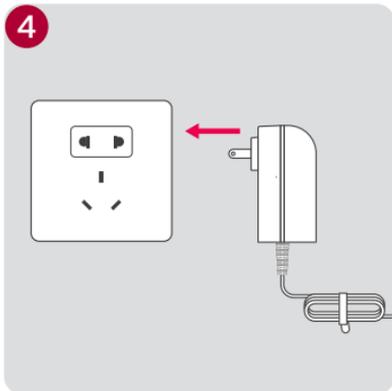

Connect the Power


- Connect the power adapter to the device.
- Plug the other end of power adapter into a power outlet.

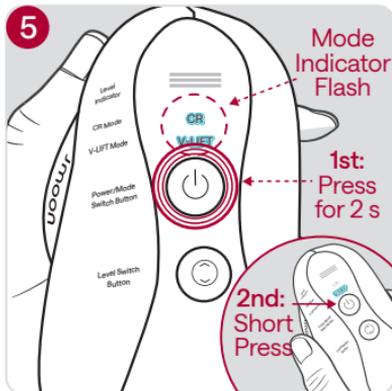

Power On (Power/Mode Switch Button)


- 1st Step:
Press and hold the Power/Mode Switch Button for 2 seconds to power on. The mode indicator will flash.
- 2nd Step:
Short press the Power/Mode Switch Button again, and the Electrode head light turn on, you can start operation from V-LIFT Mode.

* If you don't press Power/Mode Switch Button, the device will shut down automatically after 3 minutes.

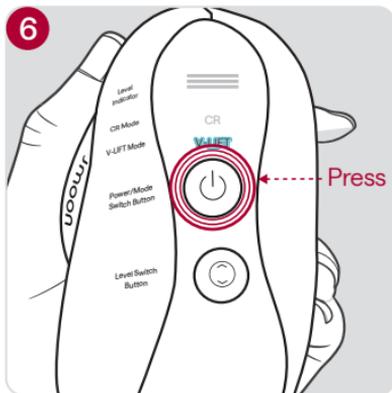


Select Mode (Power/Mode Switch Button)

Press the Power/Mode Switch button to switch modes.

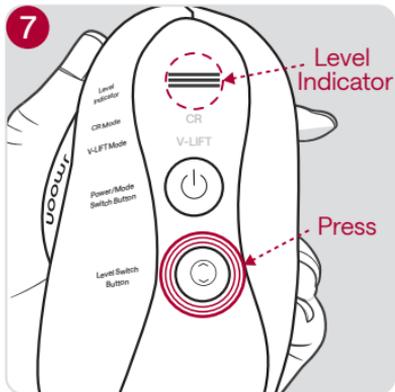

Select Level (Level button)


- Each mode has 3 intensity levels. Press Level Button to change level.
- Higher level = higher energy output.
- For 1st time use, we recommend starting from Level 1. You can adjust your level based on your comfort.

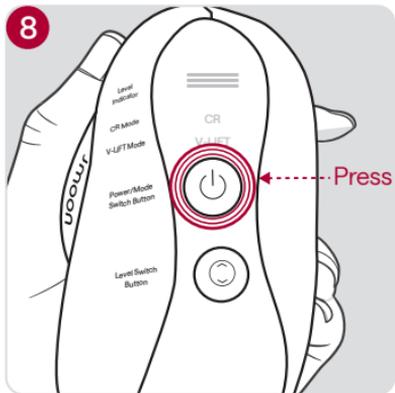

Power Off (Power/Mode Switch Button)


Press the Power/Mode Switch Button for 2 seconds to power off.

V-LIFT Mode

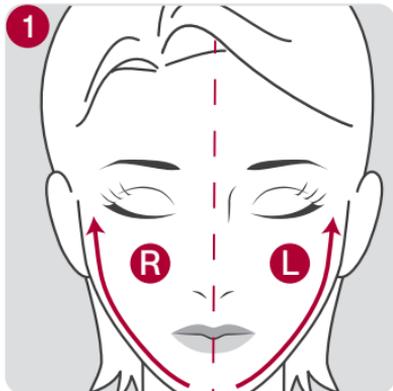


Usage Method

- 6-minute face session.
- Left and right areas: 3 minutes each.
- It's recommended to start from left face.
- Follow the arrows shown to glide along the face.

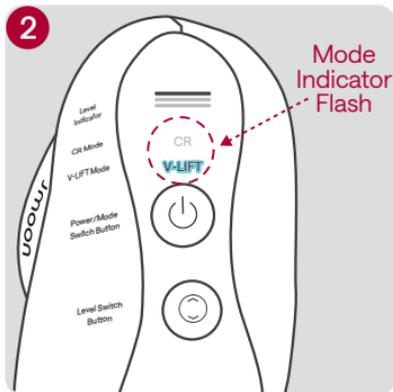

Switching Areas Reminder


- 2 beeps and 1 short vibration to remind Switch from left areas to right areas after 3 minutes
- V-LIFT Mode indicator flashes
- Press the Power/Mode Switch Button to continue operation after you hear Reminder.

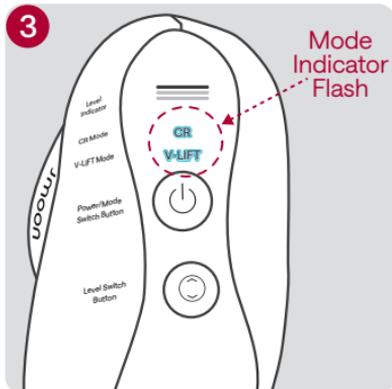

End of Session

- After the operation on right areas (3 minutes), 4 beeps and **1 long vibration indicate the end of session.**
- V-LIFT Mode and CR Mode indicator flash.
- To continue, press the Power/Mode Switch Button again.

CR Mode (Care & Renewal Mode)

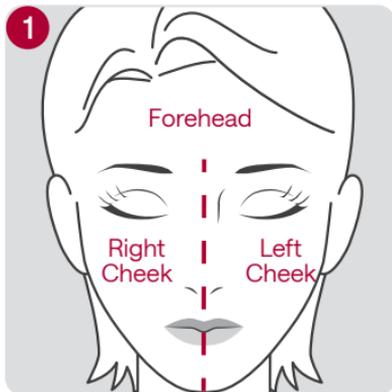


- 8-minute face session
- Recommended operation cycle: Left Cheek→Right Cheek→Forehead
- Operation time for the below 3 areas
 - Left cheek: 3 minutes
 - Right cheek: 3 minutes
 - Forehead: 2 minutes




Usage Method

- Follow the arrows shown to glide along the face.

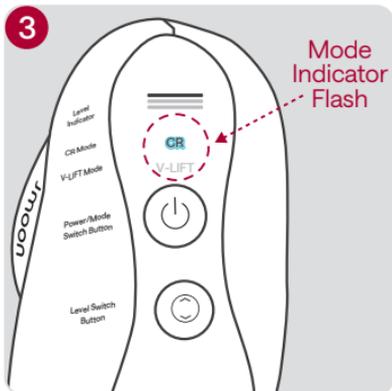

Switching Areas Reminder

- 2 beeps and 1 short vibration to remind Switch areas
- 2 times switching areas reminder
- 1st time: Remind to switch from left cheek to right cheek after 3 minutes
- 2nd time: Remind to switch from right cheek to forehead after 3 minutes
- **Press the Power/Mode Switch Button to continue operation after you hear Reminder.**

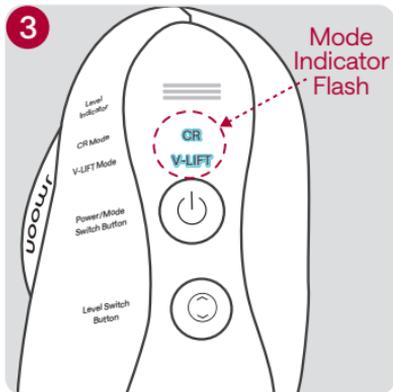

End of Session

- After the operation on forehead(2 minutes), 4 beeps and 1 long vibration indicate the end of session.
- V-LIFT Mode and CR Mode indicator flash.
- To continue, press the Power/Mode Switch Button again.

Treatment Schedule

1st day	2nd day	3rd day	4th day	5th day	6th day	7th day
V-LIFT	CR	V-LIFT	CR	V-LIFT	CR	Rest

Choose the appropriate mode based on your skin condition at different times

1. CR mode: 8 minutes per session (3 minutes on the left cheek, 3 minutes on the right cheek, and 2 minutes on the forehead); 3-4 sessions each cycle for better results.
2. V-LIFT mode: 6 minutes per session (3 minutes on the left cheek and 3 minutes on the right cheek); 3-4 sessions each cycle for better results.
3. **To avoid over-treatment and protect your skin barrier, limit each session to a maximum of 15 minutes** (excluding the time for applying gel or other preparations).

CONSISTENCY IS KEY.

Make sure you use your device consistently for 4 weeks to achieve results. Snap before-and-after photos to easily track your progress instantly and over time.

Note: Individual results may vary.

